




**BET
ACCESSORY
PRODUCTS**

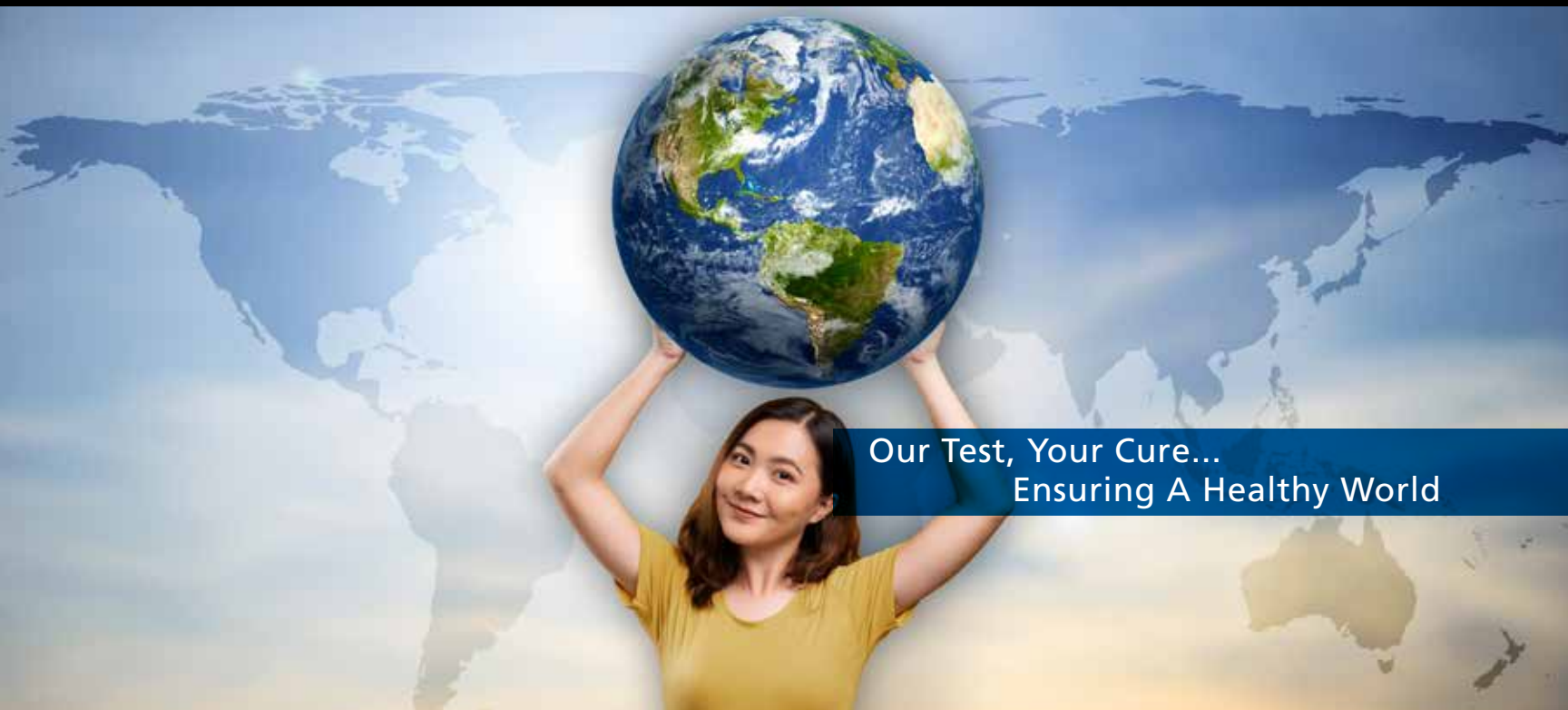

**Our Test, Your Cure...
Ensuring A Healthy World**

Your Endotoxin Experts!

Experts in Endotoxin and Glucan Detection for nearly 50 Years.

Disposable Products

Pyroclear®

Pyroclear® brand products are the first disposables in the industry that are certified to be free of interfering endotoxin and (1→3)-β-D-Glucan contamination. Pyroclear® products include depyrogenated test tubes, 96-well microplates, pipette tips and LAL Reagent water. These products are designed to reduce Out of Specification (OOS) investigations due to contaminated consumables.

LAL Reagent Water (LRW)

LRW is intended for reconstitution of BET reagents and CSE, and to dilute samples and standards for BET assays. LRW is not for human or animal injection. LRW contains less than 0.001 EU/mL endotoxin and less than 1.56 pg/mL glucan.

#W0051-10 5.5 mL/bottle, 10 bottles/pack

#W020P 20 mL/bottle, 10 bottles/pack

#WP050C 50 mL/Plastic bottle, 30 bottles/pack

#WP100C 100 mL/Plastic bottle, 30 bottles/pack

#WP500C 500 mL/Plastic bottle, 12 bottles/pack

#WP1000C 1L/Plastic bottle, 12 bottles/pack

Pyrotubes®

#TK100-10 8 x 75 mm borosilicate glass for Pyros Kinetix® Flex, 50 tubes/pack, 10 packs/carton

#TS050-10 10 x 75 mm soda lime glass for gel-clot method, 52 tubes/pack, 10 packs/carton

#TB050-5 10 x 75 mm borosilicate glass for turbidimetric method, 52 tubes/pack, 5 packs/carton

#TB240-5 12 x 75 mm borosilicate glass (for dilutions only), 42 tubes/pack, 5 packs/carton

#TB013-5 13 x 100 mm borosilicate glass, 18 tubes/pack, 5 packs

#TB16C 16 x 90 mm depyrogenation tubes with aluminum caps, 65/pack

Pyroplates®

#CA961-10 96-well microplate, 10 pack

#CA961-50 96-well microplate, 50 pack

Precision Pipette Tips

#PPT25 250 µL tips, 96 tips/box, 10 boxes/pack

#PPT10 1000 µL tips, 96 tips/box, 8 boxes/pack

Reconstitution Buffers

Pyrosol® LAL Reconstitution Buffer

Pyrosol® is an FDA licensed buffer for reconstituting Pyrotell® Multi-Test Vial or Pyrotell®-T reagents. It is used when testing electrolytes, strongly buffered solutions (especially bicarbonate buffers), and solutions for which it is difficult to adjust pH into the required range. Pyrosol® buffer is also available with a pH indicator for gel-clot applications.

#BR051-5 Pyrosol® Buffer with pH indicator (gel-clot method only) 5.5 mL/vial (5 pack)

#BR051-25 Pyrosol® Buffer with pH indicator (gel-clot method only) 5.5 mL/vial (25 pack)

#BC051-5 Pyrosol® Buffer without pH indicator 5.5 mL/vial (5 pack)

#BC051-25 Pyrosol® Buffer without pH indicator 5.5 mL/vial (25 pack)

#BC554-1 Pyrosol® Buffer without pH indicator 55 mL/vial (1 pack)

Glucashield® (1→3)-β-D-Glucan Inhibiting Buffer

Glucashield® buffer is used to reconstitute LAL and render the reagent insensitive to (1→3)-β-D-Glucan interference by effectively blocking the Factor G pathway of the endotoxin clotting cascade. For use with Pyrotell® Multi-Test Vials, Pyrotell®-T, Pyrochrome® and Chromo-LAL.

#GB051-5 Glucashield® Buffer 5.5 mL/vial (5 pack)

#GB051-25 Glucashield® Buffer 5.5 mL/vial (25 pack)

Control Standard Endotoxin (CSE)

Control Standard Endotoxin (CSE) is a standard for endotoxin testing. It is a purified extract from *E. coli* O113:H10, the same strain used for the United States Pharmacopeia and the European Pharmacopeia Reference Standard Endotoxin (RSE).

CSE is an economic alternative to the RSE. CSEs are standardized against the RSE as indicated on the Certificate of Analysis, so that results can be reported in Endotoxin Units (EU) and International Units (IU). CSE can be used for all routine BET testing. A 10 ng/vial CSE is made specifically for use with our Pyrochrome® chromogenic reagent.

Performing the Test

CSE is used to make standard curves and controls when performing the BET assay. The concentrations used are dependent on the type of assay and for photometric methods (chromogenic and turbidimetric), the detection range required.

Reconstitution

CSE is reconstituted with LAL Reagent Water (LRW). Please refer to the Certificate of Analysis when using CSE. A Certificate of Analysis for each CSE-LAL lot combination can be obtained from ACC or online at www.acciusa.com.

Depyrogenation Controls

In addition to their use as standards for controlling BET tests, the 0.5 µg and 125 µg CSEs can be used for validation of depyrogenation processes. They may be used directly, without reconstitution, as depyrogenation indicators (recommended for 0.5 µg) or can be reconstituted and endotoxin added to challenge articles (recommended for 125 µg).

Stability

Once reconstituted, CSE stored at 2–8°C is stable for maximum storage time for the different CSE preparations as listed below:

10 ng/vial	7 days
0.5 µg/vial	4 weeks
125 µg/vial	3 months

CSE should not be frozen.

Product Benefits

- CSE is a stable preparation of endotoxin that can be used in all BET testing
- CSE E0005-1 and E0125-1 can be used for depyrogenation studies
- Certificates of Analysis for each CSE-LAL lot pairing gives a potency that is specific to the unique lot combination
- CSE can be reconstituted to achieve specific endotoxin concentrations

Control Standard Endotoxin, *Escherichia coli* O113:H10

#E0005-1	0.5 µg/vial (1 pack)
#E0005-5	0.5 µg/vial (5 pack)
#E0125-1	125 µg/vial (1 pack)
#E0125-5	125 µg/vial (5 pack)
#EC010-5	10 ng/vial (5 pack)

ACC'S PYROCLEAR® BRAND
BET ACCESSORY PRODUCTS
ARE CERTIFIED TO BE FREE
OF INTERFERING ENDOTOXIN
AND (1→3)-β-D-GLUCAN
CONTAMINATION.





Ordering Information

Customer service representatives are available to assist you with orders, pricing requests and Certificates of Analysis.

Method of Payment For United States

- Check (in US dollars) made payable to Associates of Cape Cod, Inc.
- Wire Transfer (contact Accounts Receivable for routing information)
- Credit Card (AMEX®, VISA®, MasterCard®)
If payment is to be made by credit card, the following information is required: Type of Credit Card, Card Number, Credit Card Security Code, Expiration Date of Card, and Name (as it appears on the card).

Additional Information

ACC reserves the right to institute, modify or discontinue credit limits provided to customers at any time for any or no reason.

The use of credit cards for payment may incur a fee, please see our website for ACC's policy on credit card usage at www.acciusa.com/pdfs/acctc.pdf

Outside the U.S.

Please contact your local office for information regarding method of payment. For a list of your country specific distributors please visit www.acciusa.com.

Product listings, information and fill sizes are subject to change at any time without prior notice.

Corporate Headquarters

Associates of Cape Cod, Inc.

124 Bernard E. Saint Jean Drive
East Falmouth, MA 02536-4445 USA
Tel: (800) LAL-TEST (525-8378) or (508) 540-3444
Fax: (508) 540-8680
www.acciusa.com
Customer Service: custservice@acciusa.com
Technical Service: techservice@acciusa.com
Contract Test Service: testservice@acciusa.com

United Kingdom

Associates of Cape Cod Int'l., Inc.
Deacon Park, Moorgate Road
Knowsley, Merseyside, L33 7RX, United Kingdom
Tel: (44) 151-547-7444
Fax: (44) 151-547-7400
E-mail: info@acciuk.co.uk
www.acciuk.co.uk
UK Customer Service: customerservices@acciuk.co.uk
Company Registration Number: BR002906

Europe

Associates of Cape Cod Europe GmbH
Opelstrasse 14
D-64546 Mörfelden-Walldorf, Germany
Tel: (49) 61 05-96 10 0
Fax: (49) 61 05-96 10 15
E-mail: service@acciusa.de
www.acciusa.de

All Products and Services listed herein are offered exclusively under Associates of Cape Cod, Inc.'s Terms and Conditions of Sale which can be found at www.acciusa.com/pdfs/acctc.pdf



ASSOCIATES OF
CAPE COD
INCORPORATED

Specialists in Endotoxin and Glucan Detection

Corporate Headquarters

Associates of Cape Cod, Inc.

124 Bernard E. Saint Jean Drive
East Falmouth, MA 02536-4445 USA

Tel: (888) 395-ACC1(2221) or
(508) 540-3444

Fax: (508) 540-8680

www.acciusa.com

Customer Service:

custservice@acciusa.com

Technical Service:

techservice@acciusa.com

Contract Test Service:

testservice@acciusa.com

United Kingdom

Associates of Cape Cod Int'l., Inc.

Deacon Park, Moorgate Road
Knowsley, Liverpool L33 7RX

United Kingdom

Tel: (44) 151-547-7444

Fax: (44) 151-547-7400

E-mail: info@acciuk.co.uk

www.acciuk.co.uk

Company Registration Number:

BR002906

Europe

Associates of Cape Cod Europe GmbH

Opelstrasse 14

D-64546 Mörfelden-Walldorf

Germany

Tel: (49) 61 05-96 10 0

Fax: (49) 61 05-96 10 15

E-mail: service@acciusa.de

www.acciusa.de

Associates of Cape Cod, Inc. -
a Seikagaku Group Company